

ILT750 Bili Light Meter


The ultimate light measurement solution for your neonatal jaundice treatment system

Standard Features

- Dual units with broad sensitivity range: .08 to 8000 $\mu\text{W}/\text{cm}^2/\text{nm}$, 3.0 to 300,000 $\mu\text{W}/\text{cm}^2$
- 5-digit display resolution
- Peak response at 470 nm to match bilirubin action spectrum
- ISO17025 Accredited and NIST Traceable calibration
- Class 1 Medical Device
- Cosine corrected sensor for accurate spatial response
- Band Irradiance ($\mu\text{W}/\text{cm}^2/\text{nm}$) or Irradiance (W/cm^2)
- Works with traditional sources as well as LED
- Works on bili blankets, pads or overhead lamps
- Simple to use device
- Measures, calculates and provides irradiance data on a bili blanket in a matter of seconds
- Easy to clean housing
- Compact, easy-to-hold, ergonomic design
- Uses AA batteries

VERSATILE DETECTOR

The unique design of the ILT750 detector allows you to use it in any orientation. Use the meter with the detector head attached or detached, in the face up or face down orientation, always allowing you to see the screen while measuring. Detector head cable is nearly one meter long, for easy reach into incubators

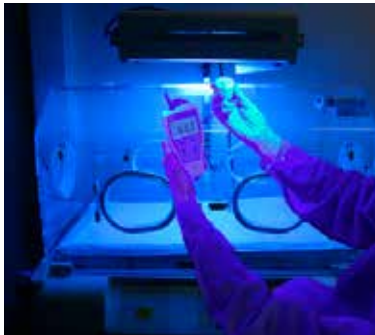

DESCRIPTION


The fastest, simplest way to validate your neonatal jaundice light source

The ILT750 completely changes the way bili light sources are measured. Say goodbye to taking multi-point measurements and doing math. The ILT750 automatically calculates readings while you seamlessly run the detector over the blanket or under the lamp. The meter will automatically give you the Min/Max and Average in a matter of seconds. No calculating or calculator required!

The ILT750 was designed with the OEM and distributor in mind. In addition to the standard features, customers can choose from an extensive list of customizable options including a custom label, start up/splash screen, custom measurement speed, digits of resolution, units of calibration, custom filtration, and more. Whether you're looking to add a light meter to your bili treatment system, or upgrading the one you have, the ILT750 will set you apart with a simple-to-use, easy-to-clean, Class 1 Medical light meter.

ILT has been a leader in light measurement technology for over 55 years. We have partnered with OEM's in many industries including medical, aerospace, industrial, manufacturing and more. The ILT750 light meter offers OEMS a fully branded, end-to-end, blue light measurement system, backed by a trusted partner in light measurement.

SPECTRAL RESPONSE

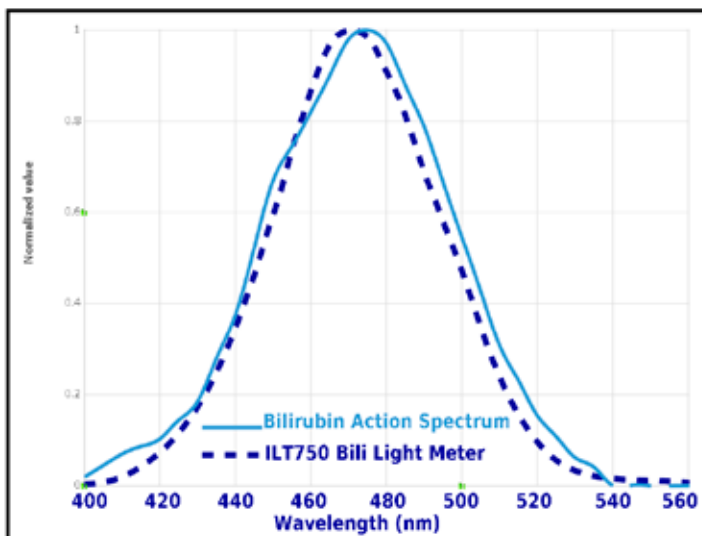

Custom filtration designs are also available to match the light source (LED, traditional lamp, or the action spectra)



InternationalLight
TECHNOLOGIES



深圳为尔康科技有限公司 联系人：曾祥满 手机：13632925349

QQ：274798107 电话：0755-28896837 地址：深圳市龙岗区沙平北路111号6008

网址：www.medicalQC.com 邮箱：szchina1718@163.com